





VELA Tango

VERSION: 600 & 600ES

Manual. nr. 105671

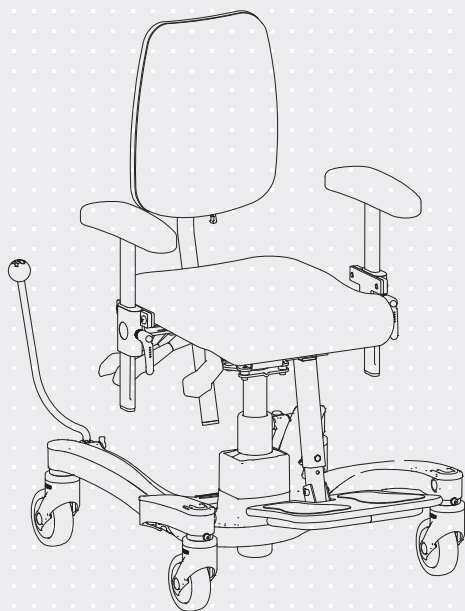

VELA Tango 600

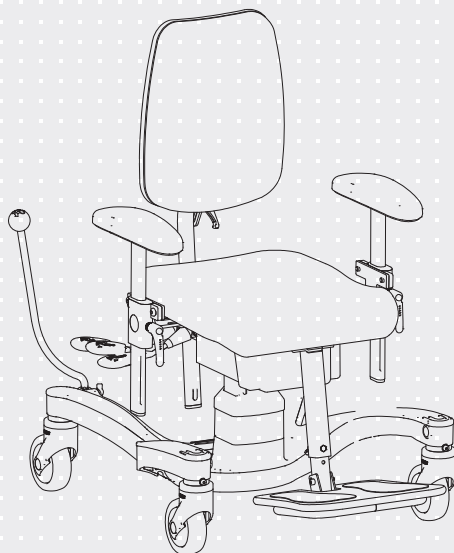

VELA Tango 600ES



ISO 9001
ISO 14001

The chair's design varies depending on the model

CHAIR FUNCTIONS



Backrest height



Backrest angle



Seat depth



Brake



Seat angle



Tilt



Armrest width




Armrest height



Seat rotation (only 600S)
(handle on left side)



Electric height adjustment

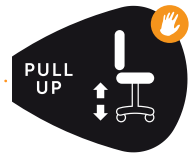

Manual height adjustment



TABLE OF CONTENTS

| | |
|--|-----------|
| 1.0. INTRODUCTION | 4 |
| 1.1. WARNINGS | 5 |
| 1.2. TROUBLESHOOTING | 6 |
| 1.3. VELA TANGO WITH ELECTRIC HEIGHT ADJUSTMENT | 7 |
| 1.3.1. Power input | 7 |
| 1.3.2. Emission | 7 |
| 1.3.3. Surge protection | 7 |
| 1.4. LABELS USED ON THE CHAIR | 8 |
| 1.5. SAFETY | 8 |
| 1.5.1. Intended use | 8 |
| 1.5.2. Transport | 9 |
| 1.6. WARRANTY | 9 |
| 1.7. UNPACKING THE CHAIR | 9 |
| 1.8. MAINTENANCE | 9 |
| 1.8.1. Wheels | 10 |
| 1.8.2. Frame | 10 |
| 1.8.3. Padding | 10 |
| 1.8.4. Recycling guide | 10 |
| 1.9. RECYCLING AND DISPOSAL | 10 |
| 2.0. SETTINGS | 11 |
| 2.0.1. Attaching the backrest | 11 |
| 2.0.2. Height adjustment (Electric) | 11 |
| 2.0.3. Soft touch contact (Electric)* | 11 |
| 2.0.4. Adjustment of seat height (Manual) | 11 |
| 2.0.5. Adjustment of seat angle | 12 |
| 2.0.6. Adjustment of seat tilt | 12 |
| 2.0.7. Adjustment of backrest angle | 12 |
| 2.0.8. Adjustment of seat depth | 12 |
| 2.0.9. Adjustment of back height | 13 |
| 2.0.10. Back latch | 13 |
| 2.0.11. Adjustment of armrest, height | 13 |
| 2.0.12. Adjustment of armrest, width | 13 |
| 2.0.13. Footrest* | 14 |
| 2.0.14. Rotation (Manual)* | 14 |
| 2.0.15. Brake | 14 |
| 2.0.16. Foot pedal for braking | 14 |
| 2.0.17. Setting the angle of the brake | 15 |
| 2.0.18. Strolling bracket* | 15 |
| 2.0.19. Charging (Electric)* | 15 |
| 2.0.20. Hand box (Electric)* | 15 |
| 2.0.21. Electric brake (Electric)* | 16 |
| 2.0.22. Folding footrest* | 16 |
| 2.0.23. 4-wheel brake* | 16 |
| 2.0.24. Accessories | 16 |

1.0. INTRODUCTION

DEAR CUSTOMER

Congratulations on your new VELA chair! We are sure this chair will give you daily pleasure. This user manual provides useful information on adjusting, operating and maintaining your chair.

Important

Please read this user manual carefully and keep it safe for when you need to refer to it. This manual is also available on VELA's website: www.vela.eu, where you can also find other relevant information.

If you have questions, experience errors or adverse events, please contact your VELA supplier.

We reserve the right to make changes.

VELA

Gøteborgvej 8-12
DK-9200 Aalborg SV
Denmark
www.vela.eu

1.1. WARNINGS



| | |
|----------------|--|
| Use | In case of an allergic reaction, stop using the chair and consult a doctor |
| Use | Only use original spare parts from VELA |
| Use | Only use the supplied charger for charging the chair |
| Use | Do not make any changes to the accessories that come with the chair |
| Use | Only use the chair for the intended purpose |
| Use | Only adjust the seat angle and backrest angle while the user is sitting in the chair and the user's feet are flat on the floor |
| Use | Always have the chair brake on when you sit down or get up |
| Use | The brake rod may only be used to operate the brake. If putting too much load on the brake rod, the brake rod/holder may brake and thereby lose its function |
| Use | On smooth surfaces such as highly polished floors or marble tiles, the chair may slide. Here, it is recommended to use a 4-wheel brake |
| Use | Never use the chair on surfaces with level differences as the chair may become unstable |
| Use | Maintenance other than cleaning should always be performed by a professional |
| Use | Never service the chair while it is in use |
| Use | Always follow the recommended maintenance and cleaning instructions |
| Use | The chair's functions are intended for the user. Only the user or trained professionals should adjust the chair while the user is sitting in it |
| Use | Never use the electric height adjustment for more than two minutes at a time |
| Use | In the event of an electrical error or unintentional height adjustment, push both of the buttons for height adjustment simultaneously for 10 seconds |
| Use | Risk of shearing or squeezing fingers in moving parts. Therefore, never use the chair if someone is touching underneath the seat |
| Choking hazard | Choking hazard on small parts. Never allow children to play with the chair, as small parts are removable |
| Choking hazard | If the chair is equipped with electric height adjustment, cables will be present. Never let children play with these |
| Electric noise | Using cables other than the standard ones may increase electrical emission or decrease electrical immunity |
| Electric noise | Avoid placing the chair close to appliances that emit electrical noise |
| Electric noise | The chair may only be used with accessories mentioned in this user manual, as electromagnetic interference may otherwise occur |
| Electric noise | If the chair moves unintentionally or does not work as expected, find another location for the chair |

1.2. TROUBLESHOOTING

| Error | Possible cause and solution |
|---|---|
| Electric height adjustment does not work | <p>The plug is not in the armrest – plug it in</p> <p>There is no power on the chair's battery. Charge the battery with the supplied charger</p> <p>The fuse is broken. See section 1.3.3</p> <p>The height adjustment needs to be reset. Press and hold the up and down buttons (simultaneously) for at least 20 seconds until the lift is in the lowest position.</p> |
| Manual height adjustment does not work | <p>There is too much load or uneven load on the chair. Try to adjust the height without sitting in the chair. If this does not work, the gas cartridge is defective and should be replaced. Contact your VELA supplier</p> |
| The brake function does not work | <p>The brake is not activated. See section 2.0.15., 2.0.16., 2.0.21. or 2.0.23.</p> |
| The armrest seems loose | <p>One of the hand screws might be loose. Tighten the hand screw</p> |
| The chair rolls badly | <p>Wheels are blocked. Clean the wheels for wire residue and the like.</p> |
| Brake rod or chassis is in the way when I move in the chair | <p>You have a chair with seat rotation. Use seat rotation to straighten the chair. See section 2.0.14.</p> |
| Brake rod seems loose | <p>The hand screw holding the brake rod might be loose. Tighten the screw. See section 2.0.17.</p> |
| The backrest arm slides down | <p>Tighten the hand screw further, as shown in section 2.0.1. You can tighten the hand screw hard without damaging the product.</p> |

If the above does not solve your problem, you should contact your VELA supplier.

1.3. VELA TANGO WITH ELECTRIC HEIGHT ADJUSTMENT

1.3.1. POWER INPUT

Voltage: 24Vdc

Current: 0.5A

Power: 12VA

1.3.2. EMISSION

A VELA chair with electrical functions has been tested and complies with current requirements regarding electromagnetic radiation (EMC requirements). The chair is tested to the following specifications:

:: Emission class CISPR 11, Group 1, Class B

See more details at www.vela.eu

1.3.3. SURGE PROTECTION

A fuse is installed on the chair which secures the chair against power surges. The replaceable part of the fuse has the following specifications:

:: Flat fuse, type: ATOF blade fuse, fast-acting

:: Voltage: 32V




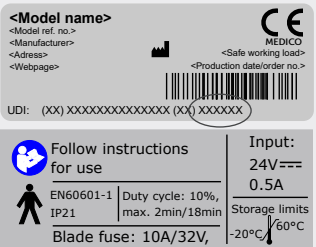
:: Current: 10A

:: Max. current: 1000A

Replacement of the fuse may only be performed by professionals, as the cause of the breach must be documented for safety reasons.

IN GENERAL

1.4. LABELS USED ON THE CHAIR

| | | |
|---------|--|---|
| Warning | Risk of shear and squeeze: Keep hands and fingers away from moving parts |  |
| Icon | Read "Instruction For Use" before use |  |
| Icon | Part used: Type B cf. IEC60601-1 covering the following areas of the chair: |  |
| Labels | 60 x 22,5 mm The last six digits are the product ID number: This label only applies to chairs with electrical functions: The chair is IP21 approved, which corresponds to protection against water droplets, as well as against touching dangerous parts with a finger. |  <p> <Model name> <Model ref. no.> <Manufacturer> <Address> <Webpage> </p> <p> CE MEDICO <Safe working load> <Production date/order no.> </p> <p> UDI: (xx) xxxxxxxxxxxxxxxx (xx) xxxxxx </p> <p> Follow instructions for use EN60601-1 IP21 Duty cycle: 10%, max. 2min/18min Blade fuse: 10A/32V, </p> <p> Input: 24V === 0.5A Storage limits -20°C / 60°C </p> |

Labels on the chair should be read at a distance of approx. 40 cm under normal daylight conditions.
For labels related to the operation of the chair, please refer to section 2.0.

1.5. SAFETY

This chair is CE marked and complies with applicable EU regulations. The chair is manufactured by: VELA, Gøteborgvej 8-12, 9200 Aalborg SV, Denmark.

The VELA Tango with electrical functions meets all safety requirements according to IEC 60601-1, including requirements for EMC according to IEC60601-1-2.

1.5.1. INTENDED USE

This device is intended to support a child in a sitting posture while performing indoor activities of daily living, including school activities.

- :: The chair is not built to stand on
- :: The wheel brake must always be activated when sitting down or standing up

All variants of this chair are tested and approved for a maximum load of 80 kg. For optimal use, the chair's recommended user weight is 15-50 kg. **Note:** Persons using a harness in the chair must always be supervised when the harness is in use.

1.5.2. TRANSPORT

When transporting the chair by car, airplane or any other means of transport, always use the chair brake. The chair is not approved for use as a passenger seat in any means of transportation and may not be used as such.

The chair is designed to be transported under the following environmental conditions:

Temperature: - 20°C to + 60°C

Humidity: 15-95%

Allow the chair to stand for 24 hours before use if it has been stored below - 10°C.

1.6. WARRANTY

VELA grants a right of claim under applicable law in the country where the chair is acquired. The right of complaint applies only when using original spare parts and accessories, as well as customisations made by professionals.

Note: if an exclusive contractual agreement has been made directly with VELA, this applies!

VELA assumes no liability for damages to the product or the user caused by:

- :: Transportation
- :: Misuse
- :: Normal wear and tear
- :: Wrong use
- :: Use of spare parts and accessories that are not original VELA products
- :: Adjustments made by unauthorised service technicians

1.7. UNPACKING THE CHAIR

The chair comes in environmentally friendly and recyclable packaging.

In case of defects or transport damage to the chair, please contact your VELA supplier right away. In such cases, do not attempt to use the chair or to fix it yourself.

On delivery you have received:

- :: VELA chair
- :: Instructions for use
- :: Charger (only for chairs with electrical functions)

Note: see item 2.0.1 for correct attachment of backrest.

1.8. MAINTENANCE

We recommend an annual inspection of the chair to identify any defective or worn parts. All removable parts should be tightened regularly. You can find various guides on www.vela.eu.

1.8.1. WHEELS

We recommend you regularly check the wheels and remove any residual threads, etc. Certain floor types accumulate dirt on the castors. Therefore, clean the castors as needed.

1.8.2. FRAME

Clean with lukewarm water and if necessary, add some pH-neutral detergent.

1.8.3. PADDING

Fabric: Clean by chemical cleaning, with furniture cleaner or with water using a hard-wring cloth.

Synthetic fabric or leather: Clean with water or mild soap using a hard-wring cloth.

Leather: Clean by vacuuming with a soft brush. Avoid using detergent.

Note: On www.vela.eu you can find cleaning instructions for all our standard upholstery types.

1.8.4. RECYCLING GUIDE

Your VELA chair can be reused by another person when you no longer need it as long as the chair is serviced by a technician beforehand.


The chair's life expectancy is five years for regular and proper use. We recommend that VELA, your VELA supplier or a trained professional perform an annual inspection with the replacement of worn parts.

1.9. RECYCLING AND DISPOSAL

This product contains recyclable materials. Therefore, the product must be disposed of in accordance with local regulations and not together with ordinary household waste. Proper disposal and recycling will help prevent negative effects on the environment and people.

Note: Electrical and electronic equipment and accompanying batteries contain materials, components and substances that can be harmful to human health and the environment if the waste is not handled properly.


Electrical and electronic equipment and batteries are marked with the icon shown below. This symbolises that electrical and electronic equipment and batteries must not be disposed of together with unsorted household waste but must be collected separately.

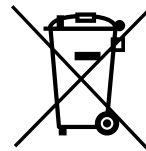

It is important that you return your used batteries to the established schemes. In this way, you help ensure that batteries are recycled in accordance with current legislation and do not cause unnecessary harm to the environment.

2.0. SETTINGS

The gas spring and the electric lifting mechanism should only be replaced by a VELA technician or a person trained in this field.

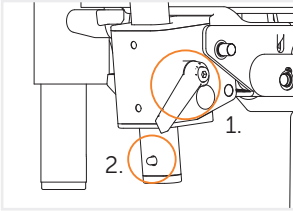

2.0.1. ATTACHING THE BACKREST

Press the locking pin spring (1) and push the backrest arm through the hole. Tighten the quick release handle (2) afterwards.

Note: The backrest arm must be pushed through completely so the pin spring locks the backrest arm. Make sure that the backrest arm cannot be pulled out again.



2.0.2. HEIGHT ADJUSTMENT (ELECTRIC)

Is operated via two switches (up/down) located on the underside of either the right or left armrest. The armrest can easily be moved to the opposite side.

Note: Height function should only be used continuously for one minute, then a nine-minute break should follow.



2.0.3. SOFT TOUCH CONTACT (ELECTRIC)*

Is operated by pressing/pulling the switch that is located at the underside of either the right or left armrest. Let go of the switch when you reached the wanted sitting height.

Note: The height function should only be used continuously for one minute, then a nine-minute break should follow.




2.0.4. ADJUSTMENT OF SEAT HEIGHT (MANUAL)

Adjust the seat height by pulling up the front lever placed at your right hand.


* Depending on model

Note: Setting the angle of backrest and seat should always be done with your feet flat on the floor.

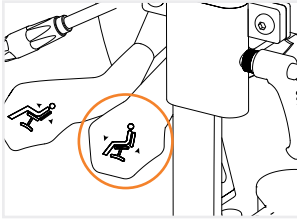

2.0.5. ADJUSTMENT OF SEAT ANGLE

Adjust the seat angle by pulling up the middle lever placed at your right hand.

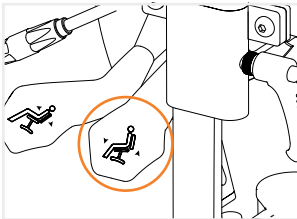

2.0.6. ADJUSTMENT OF SEAT TILT

The entire seating unit moves (seat and backrest together) by lifting the middle handle on the right side.



2.0.7. ADJUSTMENT OF BACKREST ANGLE

Adjust the back tilt by pulling up the rear lever placed at your right hand.


Only operate the handle when someone is sitting in the chair - as the backrest tilts forward.



2.0.8. ADJUSTMENT OF SEAT DEPTH


Adjusted by pulling up the rear handle on the right side (backrest angle). Adjusted so that the seat depth suits the user of the chair.

Operate the handle only when someone is sitting in the chair, as the backrest tilts forward.

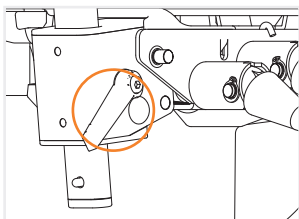


2.0.9. ADJUSTMENT OF BACK HEIGHT

Adjust by loosening the quick release handle on the backrest and regulate the height. Remember to tighten the quick release handle again.

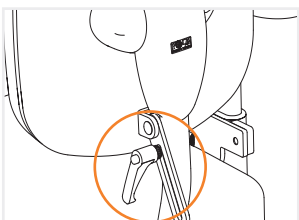


2.0.10. BACK LATCH

The backrest can be set to a specific angle by adjusting and secure the bolt which is mounted on the back of the chair back.

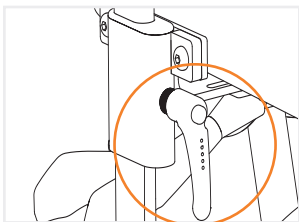

2.0.11. ADJUSTMENT OF ARMREST, HEIGHT


Adjust by loosening the screw/lever on the armrest bar and regulate the height of the armrest. Remember to tighten the screw/lever again.

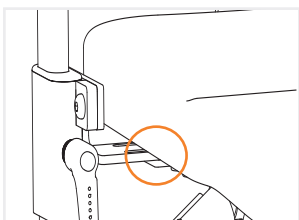

2.0.12. ADJUSTMENT OF ARMREST, WIDTH

Adjust by loosening the screw underneath the seat and regulate the width of the armrest. Remember to tighten the screw again.

* Depending on model

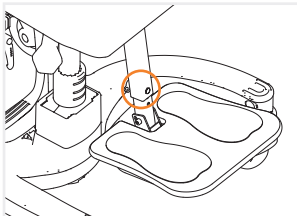

2.0.13. FOOTREST*

The Footrest can be adjusted in length by loosening the hand screw in the center of the rod on the footrest. Remember to tighten the screw afterwards. The foot plate can be folded up so that it does not interfere when getting in and out of the chair.


Warning: This footrest is not made to stand on - if this is done, there is a risk of the chair tipping over.



2.0.14. ROTATION (MANUAL)*

The seat is rotated by pulling the rear handle on the left side.


Note: Cannot be combined with Footrests, section: 2.0.13.

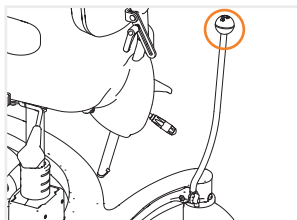

2.0.15. BRAKE


The brake is activated by moving the handle on the side of the chair forward. Unlock the chair by moving the handle back into the upright position.

Note: See Setting the angle of the brake, section 2.0.17. for adjusting the position of the brake and the distance to the seat.

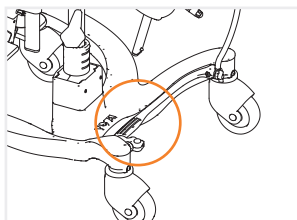

2.0.16. FOOT PEDAL FOR BRAKING


The brake can be activated with the foot pedal by pressing the foot against the brake shaft in the area shown.

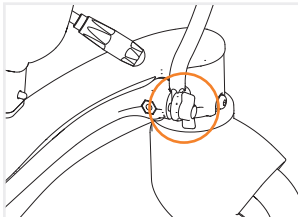

2.0.17. SETTING THE ANGLE OF THE BRAKE


Loosen the thumbscrew and remove the brake rod from the holder. Insert the brake rod in the desired side and angle. Remember to retighten the thumbscrew.

Note: The brake can be set at 30° intervals. The brake rod can be moved forwards/backwards and closer to/further from the seat as required. The brake can also be moved to the opposite side.

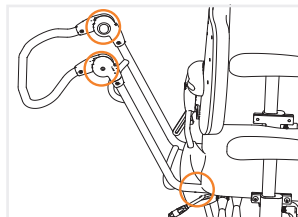

2.0.18. STROLLING BRACKET*

Strolling bracket for transporting a person in the chair. The strolling bracket is mounted to fittings under the seat using two finger screws. This is also where the depth is adjusted. The height is adjusted by pressing the two buttons on the side of the bracket while the handle is positioned at the desired height.

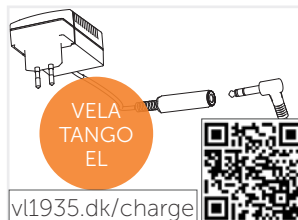

2.0.19. CHARGING (ELECTRIC)*

Charge every day Remove the operation plug from the back of the armrest and place it in the charger as shown. The charger is plugged into an electrical outlet. When the chair is charging, the LED on the charger is yellow. When the chair is fully charged, the LED turns green.

Note: The chair must not be used during charging.



2.0.20. HAND BOX (ELECTRIC)*

The hand box is used via push buttons on the front. VELA hand boxes can have one or more functions. The icons on the hand box show the button functions.

Note: Insert plug in hand box to operate functions.
Warning: The chair's operating specifications for temperature change when using the hand box to + 5°C to + 40°C!




* Depending on model



2.0.21. ELECTRIC BRAKE (ELECTRIC)*

The electric brake is activated and deactivated using the hand box or switches. It is advisable to activate or deactivate the brake completely when using this feature. Can also be used together with the 4-wheel brake (see section 2.0.23.)


Note: The brake can only be used via switches under the armrest, the hand box, or Soft touch.

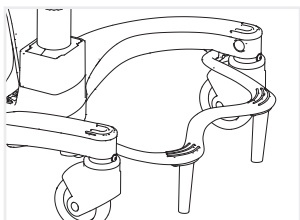

2.0.22. FOLDING FOOTREST*

Stabilising foot support that can also be used as a footboard when moving in and out of the chair.

Note: Cannot be combined with Footrest, section: 2.0.13.

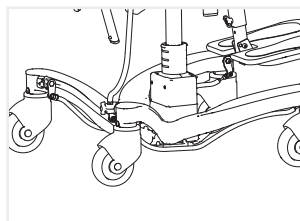

2.0.23. 4-WHEEL BRAKE*

Is activated via the hand/foot brake on all chairs except on chairs with electric brake, where it is activated via hand box or switches. 4-wheel brake brakes on all wheels and provides a stable seat on uneven/skewed surfaces.

Note: On the electric version it is advisable to activate or deactivate the brake completely when using this feature.

2.0.24. ACCESSORIES

This chair is available with a wide range of accessories. See a complete list of accessories at: www.vela.eu

Note: Charger is included as an accessory when purchasing a VELA chair with electric functions.

* *Depending on model*